

Anesthetizing Box

Cat. No. 7900 (rodents) 7910 (rabbits)

General

Our Induction Boxes are conveniently dimensioned induction boxes, featuring a sliding lid. They are made of Perspex and prove to be particularly useful to confine laboratory animals during anesthetizing.

The **7900**, for small rodents, is **dimensioned 25x13x13(h)cm**; the larger model 7910, for rabbits is **dimensioned 40x22x21(h)cm**;

The transparent acrylics permits the animal to be kept under constant observation.

Two tubing connectors of nickel plated brass are fitted into each end, one located at the top of the box and the other at the bottom.


Any (non-explosive!) gas mixture can be used. In case small quantities of liquid, as ether or chloroform are used, soak a cotton wool flock and place it in a small Becker, in-side the box.

For more demanding application, and higher safety, an airtight model, with latch, is also available, see picture below.



**Our Induction chambers
are ideal to work with
our new
Anesthesia Systems**

**TO CONFINE SMALL
LABORATORY ANIMALS
DURING
ANESTHETIZING**

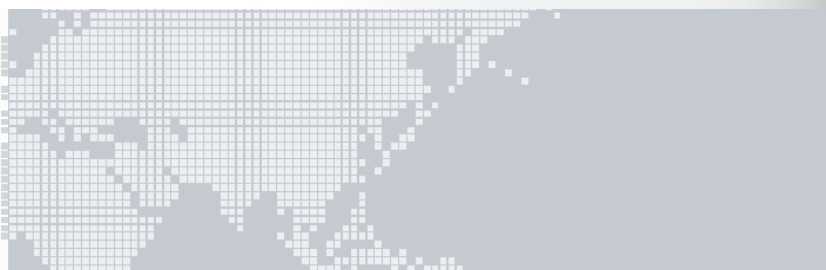

Ordering Information

- **7900** Induction box for small rodents
25x13x13(h)cm, ID 23x12x12(h) cm
- **7910** Induction box for rabbits
40x22x21(h)cm, ID 38x20x19(h) cm
- **2100-790** Airtight model, with latch, see picture
25x13x13 (h) cm, ID 21x11x13(h) cm

